


QUALITY PILL SPLITTERS AND PILL PULVERIZERS

These easy tools are ideal for elderly people, children and hospital attendants. Great also for medicating pets. They avoid dispersion of medicine and guarantee integrity of the dose. Easy to clean, also in dishwasher. Made in non-toxic, unbreakable material.

- 25729 PILL SPLITTER - bulk
- 25730 PILL SPLITTER - with box
- 25743 ULTRA PILL SPLITTER - blister

Accurately divide vitamin tablets, large and small pills for reduced dosage or easier ingestion. Simple to operate.

Top quality products, made in Europe (25729-30), USA (25743)

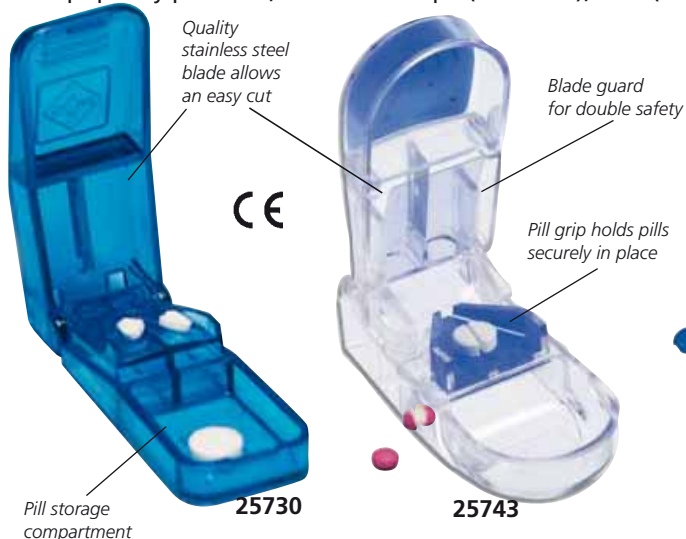


- 25747 DISPLAY FOR SPLITTER (25730) PULVERIZER (25731) - empty (english language)

Suitable for 4 pill splitters and 4 pill pulverizers. Available in any language on request.

- 25741 PILL PULVERIZER - bulk
- 25731 PILL PULVERIZER - with box
- 25744 ULTRA PILL PULVERIZER - blister

Powerful pulverizing point powders pills of any shape and size, with an easy twist of the wrist.

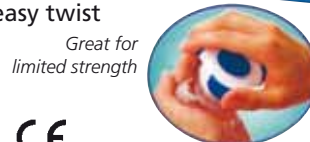

Great for limited strength



MIXED COLOUR PILL SPLITTERS

- 25748 PILL SPLITTER - box of 12 pcs
This box includes 12 pill splitters in 4 colours: light blue, green, transparent lavender and transparente orange.
- 25749 PILL SPLITTER - display box of 12 pcs
Counter display box including 12 pill splitters, each one in a single carton box in 4 colours: light blue, green, transparent lavender and transparente orange.




3 SIDES DISPLAY: FOR DIAGNOSTIC TESTS, EASY LIFE, PILLS DEVICES

METAL DISPLAY

Metal turnable display stand with 36 hangs and 3 baskets for many products. It has 3 attractive sides: for diagnostic tests (see page 4 to 7 and 107-108), for aids for daily life (see page 47) and for pill devices (see page 46 and above). Offered empty or full with a selection of 39 different GIMA products, suitable for pharmacy or shop.

Head of display language	METAL DISPLAY	
	Empty	With products
English	25520	25530
Italian	25521	25531
French	25522	25532
Spanish	25523	25533
Other on request	25529	25539

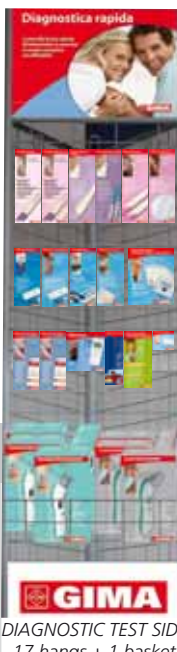


DIAGNOSTIC TEST SIDE
17 hangs + 1 basket

3 SIDE DISPLAY WITH PRODUCTS INCLUDES:

- **PILL DEVICES SIDE - 10 hangs, 1 basket**
 - 6 pill splitter (25730)
 - 6 pill pulverizer (25731)
 - 3 pill box timer (25765)
 - 6 pocket pill box (25750)
 - 6 daily pill box (25752)
 - 6 twice a day planner (25759)
 - 6 7 days planner (25760)
 - 6 burnfree® kit (34884)
 - 2 automatic sphygmomanometer (32787)
 - 2 wrist sphygmomanometer (32786)
- **EASY LIFE SIDE - 9 hangs, 1 basket**
 - 3 cutlery set - 4 pieces (28920)
 - 6 spoon & fork holder (28925)
 - 6 pencil grip - box of 2 (28901)
 - 6 key turner - set of 2 (28940)
 - 6 easy opener (28950)
 - 6 jar opener & closer (28951)
 - 6 zipper & button puller (28911)
 - 6 shopping bag carrier (28945)
 - 6 magnifying glass (28906)
 - 2 family nebulizer (28106)
- **DIAGNOSTIC TESTS SIDE - 17 hangs, 1 basket**
 - 6 pregnancy test - box of 1 (29104)
 - 6 pregnancy test - box of 2 (29107)
 - 6 menopause test - box of 2 (29109)
 - 6 ovulation test - box of 5 (29112)
 - 6 ovulation test - box of 6 strips (29114)
 - 4 breast pad (29120)
 - 4 10 parameter urine strips - 9 strips (24070)
 - 4 cholesterol strips - box of 2 (24521)
 - 4 PSA, prostate test - box of 1 (24532)

- 4 colon disease test - box of 5 (24531)
- 4 helicobacter pylori test - box of 1 (24529)
- 3 multidrug test - 6 parameters (24560)
- 3 pocket alcohol tester (24499)
- 6 digital thermometer (25560)
- 6 flexi digital thermometer (25563)
- 6 frog digital thermometer (25566)
- 4 nipple-like digital thermometer (25568)
- 2 IR ear thermometer (25580)
- 2 forehead thermometer (25585)

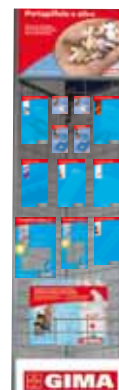

PILLS DEVICES SIDE
10 hangs + 1 basket



EASY LIFE SIDE
9 hangs + 1 basket

VIALS OPENER

• **25739 VIALS OPENER - box of 20**

A practical tool to help opening of vials in a safe way.




25739

PILL BOX AND PILL BOX TIMER

• **25732 - 7 DAYS PILL BOX**
blue translucent

Translucent top enables the user to see if medication has been taken

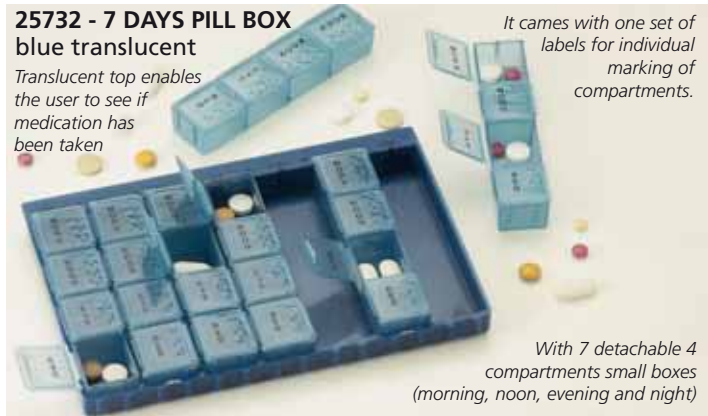

It comes with one set of labels for individual marking of compartments.

With 7 detachable 4 compartments small boxes (morning, noon, evening and night)

• **25765 PILL BOX TIMER**

Pocket pill box timer in ABS material, size 90x50x19 mm, supplied in carton box with batteries.

Double pill box space to store pills and tablets. It remains when you have to take your medication with a countdown up to 19 hours and 59 minutes.



25765



25757 - TWICE A DAY PLANNER - Italian
25759 - TWICE A DAY PLANNER - English
2 colours translucent



Available from stock with English or Italian printings.
We can print in any language, minimum 1000 pcs.



25758 - 7 DAYS PLANNER - Italian
25760 - 7 DAYS PLANNER - English
4 colours translucent



Contoured design for easy pill removal

Available from stock with English or Italian printings.
We can print in any language, minimum 1000 pcs.



Holds up to 50 pills

L M M G V S D
T T W J F S D

Supplied with a set of labels for each day suitable for 6 languages: English, French, German, Italian, Spanish, Portuguese



25750




25735



25754



25734

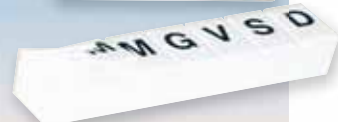

25756

- 25750 PILL BOX pack of 2 pcs - blue/pink - blister
- 25735 DAILY PILL BOX - yellow - box
- 25752 DAILY PILL BOX - blue translucent - blister

- 25754 DAILY PILL BOX - white - blister
- 25734 WEEKLY PILL BOX - white - box
- 25756 WEEKLY PILL BOX - white - blister

	25750	25735/25752	25734	25732	25757/25759	25758/25760	25754	25756
English language printed on pill boxes	yes	-	-	yes	yes	yes	-	no
Labels in other languages*	-	-	yes	yes	yes	yes	-	-
Colour	blue/pink	yellow/blue	white	translucent	translucent	translucent	white	white
Compartments	2	6	7	28(4x7)	14(2x7)	28(4x7)	4	7
Size cm	7x5x2	8x5.5x1.5	13x2.5x2	17x12x2	18.5x9.5x3	21x14.5x3.5	95x2x2	15,5x3x2
Package	blister	box/blister	box	box	box	box	blister	blister

*Supplied with a set of labels with day in 6 languages (English, French, German, Italian, Portuguese, Spanish) and 1 white label for compartment for individual marking.

EASY LIFE: AIDS FOR DAILY LIVING



• **28920 CUTLERY SET**
Set of 4 pieces (knife, fork, small and large spoon). Made in stainless steel with thick, ribbed rubber handle (10 cm) for easy grip



• **28925 SPOON & FORK HOLDERS - pack of 2 pcs.**
These portable, stylish holders fit onto virtually any style spoon and fork providing a secure grip (spoon and fork are not included).



• **28901 PENCIL GRIP - pack of 2 pcs.**
Rubber grip for pens or pencils. Helpful aids for elderly or disabled people. These cushions are designed in both shape and density to alleviate stress while providing a secure grip.



• **28945 SHOPPING BAG CARRIER**
This product provides an effort-free handle that remains constantly supportive, regardless of bag loads.



• **28950 EASY OPENER**
This all in one product easily open bottle caps, pull tabs, and even those impossible box tops! Designed to provide safety for the user.



• **28951 JAR OPENER & CLOSER**
This uniquely designed product concentrates hand power to remove lids and turn facets effortlessly.



• **28915 LONG BATH SPONGE**
These lightweight sponges feature soft polyfoam attached to a plastic handle. Designed as bath sponges, Available with contoured sponge.



• **28906 MAGNIFYING GLASS**
Magnifies two times larger than original size. Large enough to do what's needed, yet small enough to be totally portable, fitting easily into purse or pocket.



• **28912 SOCK AID NOTCH**
Soft terry cloth on one side grips the sock, while smooth nylon on the other side reduces friction, allowing your foot to slip in easily. Two long loop handles are easy to hold and assist with pulling the sock onto your foot. Ideal for those who have the use of only one hand or with decreased strength or coordination.



• **28911 ZIPPER AND BUTTON PULLER**
For those with weak or arthritic fingers, there's probably no more difficult task than fastening a button! This unique product eliminates the frustration by pulling the button or zipper into position. Folds to half-length to fit easily into pocket or purse.



• **28940 EASY KEY TURNER**
pack of 2 pcs
These enlarged key turners are designed with the turning effort already begun, so completion of the turn is accomplished with minimal effort. Two different colours per package.



• **28941 MULTIPLE KEY TURNER**
The key turner holder is a heavy-duty holder designed to assist person with limited wrist strength or movement. People with arthritis or other limited hand mobility can manipulate keys with ease. House keys and most car keys (without plastic head covers) fit into the slots and are held by a single screw. Can hold 1 to 3 keys. Keys fold into handle when not in use to slip into pocket or purse.



• **28963 GRASPING TOOL - foldable**
Your extended reach grasping tool will help you to eliminate bending or climbing for those hard-to-reach items. Length: 80 cm, foldable in 2. Locking tab holds the item in place.



• **28960 GRASPING TOOL - 3 WAYS - 66 cm**
Aluminium body, PP handle. Three grasping ways: steel hook, magnetic tip, rubber head.


ANTI-DECUBITUS CUSHIONS

MADE IN ITALY




28520

PREVENTION FOR LOW DECUBITUS RISK

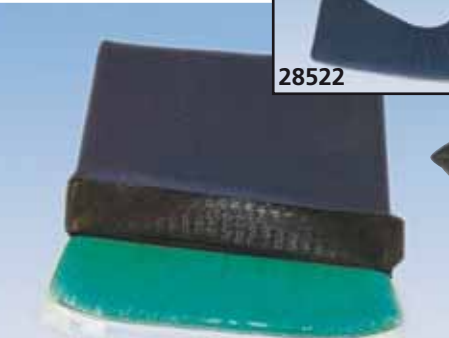

28524

PREVENTION FOR MEDIUM DECUBITUS RISK

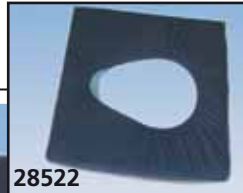

28522



28526

PREVENTION FOR MEDIUM-HIGH DECUBITUS RISK


- **28520 ANTIDECUBITUS GEL CUSHION 40x40x4 cm**
Fluid gel grants a perfect distribution of patient's weight by sensibly decreasing the pressure on capillars and avoiding their eventual occlusion. It also gives a pleasant sense of freshness.
- **28522 VISCOELASTIC CUSHION with hole 40x40x3 cm**
- **28524 VISCOELASTIC CUSHION without hole 40x40x3 cm**
It is composed by a Polyurethane gel layer with high pressure absorbing capacity and by a layer of Polyurethane expanded foam. The two layers are wrapped and protected by a special cover easily removable, with transpirant sitting side and water-proof and anti-slip bottom side. Thanks to the Polyurethane gel the contact pressure decreases and improves the distribution of the same. Cover is made of polyamide threaded with PU,

water-proof, transpirant to steam, elastic hypoallergenic and washable in washing machine up to 40°.

- **28526 AIR FLOATING CUSHION 40x40x6.5 cm**
- **28528 AIR FLOATING CUSHION 44x44x6.5 cm**
Latex cushion with communicating cells. The particular shape of the air cells allows a perfect body positioning and a correct pressure distribution. The alveolar cells allow an optimal adaptation to human shape increasing perspiration and air circulation thus reducing skin maceration. Particularly suitable for patient with skin wounds, medullary lesions, multiple sclerosis, bifid spine. The device is supplied with an anti-slip cover, a pump for inflation and a reparation kit. Available with single valve.

ELECTRIC HEATING PAD

- **28655 ELECTRIC HEATING PAD** Suitable for many applications.
- **28656 ELECTRIC "SAND" HEATING PAD**
The use of purified marine sands make it possible to naturally transmit their beneficial warmth. External cover in cotton (94%) and nylon, internal bag in PVC. Operating power: 230V-50/60 Hz - 75W (28656) - 55W (28655)

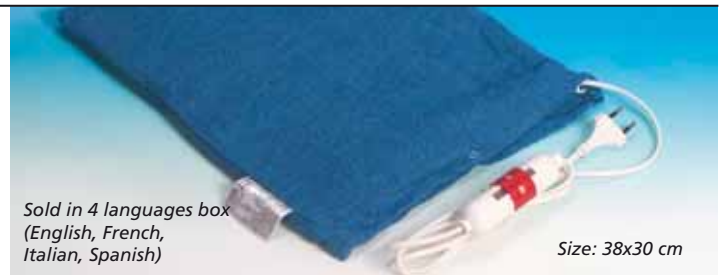

Sold in 4 languages box
(English, French,
Italian, Spanish)

Size: 38x30 cm

SOFT AND RIGID COLLARS

SOFT CERVICAL COLLARS

This collar is made of soft expanded polyurethane lined with an anti - allergic 100% cotton fabric. Velcro fixing on the back. Available in 3 sizes.

RIGID CERVICAL COLLARS

Made of polyethylene, with comfortable resting support in vinyl leather and height adjustable. Velcro straps on the back. Available in 3 sizes.

Soft

Rigid



ARM SLINGS



Suitable for supporting a sprained or broken arm during recuperation. The shoulder strap is adjustable to allow a correct positioning of the arm. Stitched from cotton fabric requiring normal laundry care. Fits any arm. Available in 3 sizes.

GIMA code	POUCH ARM SLINGS
28620	Arm sling - small
28621	Arm sling - medium
28622	Arm sling - large

INVALID RINGS



Quality invalid rings made of latex-free rubber.

GIMA code	INVALID RINGS LATEX - FREE RUBBER
28610	Invalid ring Ø 35 cm
28611	Invalid ring Ø 40 cm
28612	Invalid ring Ø 45 cm

GIMA code	SOFT AND RIGID NECK COLLARS	CE
34600	Rigid cervical collar 32-40 cm (neck circumference)	
34601	Rigid cervical collar 41-46 cm (neck circumference)	
34602	Rigid cervical collar 47-52 cm (neck circumference)	
34610	Soft cervical collar 43 x 8 cm - small (h 7 cm)	
34612	Soft cervical collar 46 x 8 cm - medium (h 8.5 cm)	
34614	Soft cervical collar 49 x 10 cm - large (h 10.5 cm)	

RAISED TOILET SEATS




• **27743 RAISED TOILET SEAT** height 10 cm with central fixing system and armrest. Max load 130 kg.

• **27741 RAISED TOILET SEAT** height 10 cm with lateral fixing system. Max load 200 kg.

• **27739 SOFT RAISED TOILET SEAT** - height 10 cm Made in polyurethane. Size 370x105x130 mm. Max load 200 kg.

TUB RAIL

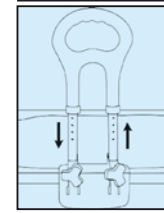

Adjust the handle by C-clips to proper height for better grip


Height adjustable 48.5-56.8 cm

• 27734 TUB RAIL

- non slip comfort grip
- easy tool-free assembly
- 2 separate clamps hold firmly to bathtubs with sloped inner walls
- handle with 5 position height adjustment
- oversized handle provides additional stability.




SHOWER AND BATH CHAIRS

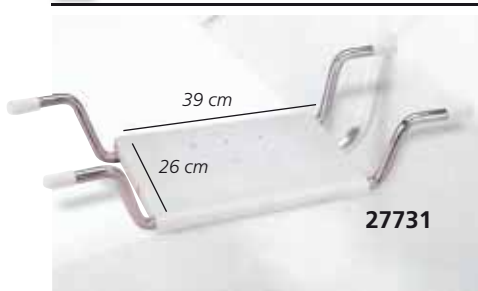

• 27731 BATH CHAIR

- clamps onto the edge of the bath, no drilling or crewing required
- swivel seat with 90 degree position lock
- handle mounts either right or left side
- suitable for most types of bath with clamp adjustable from 50 to 150 mm
- dual clamps have non-marring pads and swivels to fit shaped baths
- safety handle provides support and balance when getting in or out of bath
- non-slip seat surface/ fixed backrest
- maximum load 115 kg (250 lbs)

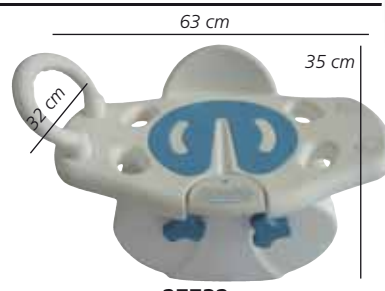

27732



27732



• 27730 SHOWER CHAIR with backrest

Seat and backrest in plastic material Height adjustable (40-50 cm), non slip end.

• 27733 ADJUSTABLE SWIVEL STOOL

- easy tool-free installation
- swivel seat with 90 degree position lock at every right angle location
- seat with 8 height adjustments, by 2.5 cm (1") increment
- non-slip central seat pad and soft under-covers
- 135 kg (300 lbs) maximum load

SUCTION GRAB BARS




• **27735 GRAB BAR** - 445 mm

• **27736 GRAB BAR** - 547 mm

• **27737 DOUBLE GRAB BAR** - 890 mm

Innovative line of GRAB BAR, can be quickly installed every where without need of drilling and screwing.

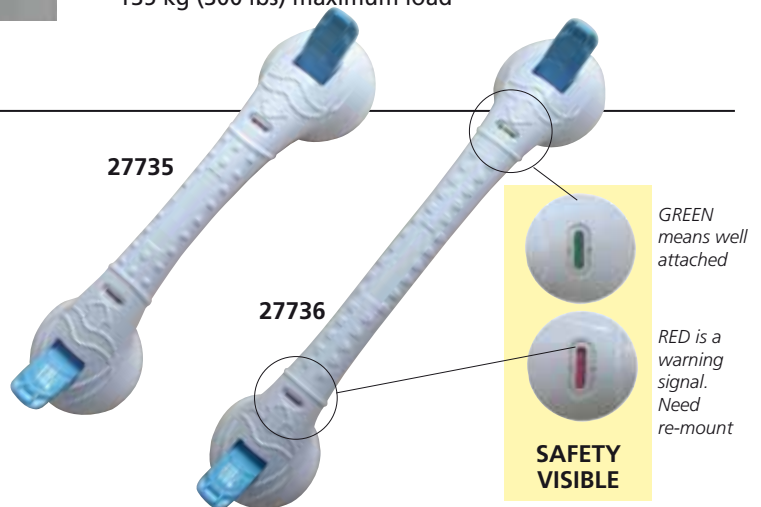

GREEN means well attached

RED is a warning signal. Need re-mount

SAFETY VISIBLE



Patented suction cups (ø 93 mm). Holds up to 60 Kg (150 lbs) lateral and 45 or 85 kg (27737) vertical weight. Made of heavy duty ABS, ideal solution to support elder people during daily bathroom use.


ROLLATORS

• 27724 ROLLATOR WITH SEAT

With seat, large front basket and brakes on wheels.
Size 71x64 cm, adjustable. Height 81-93 cm. Seat at 55 cm.
Weight 8 kg. Load: 100 Kg.

• 27721 ROLLATOR

Size 58x54 with adjustable height 51-82 cm.
Weight: 4.5 kg. Load: 100 kg.

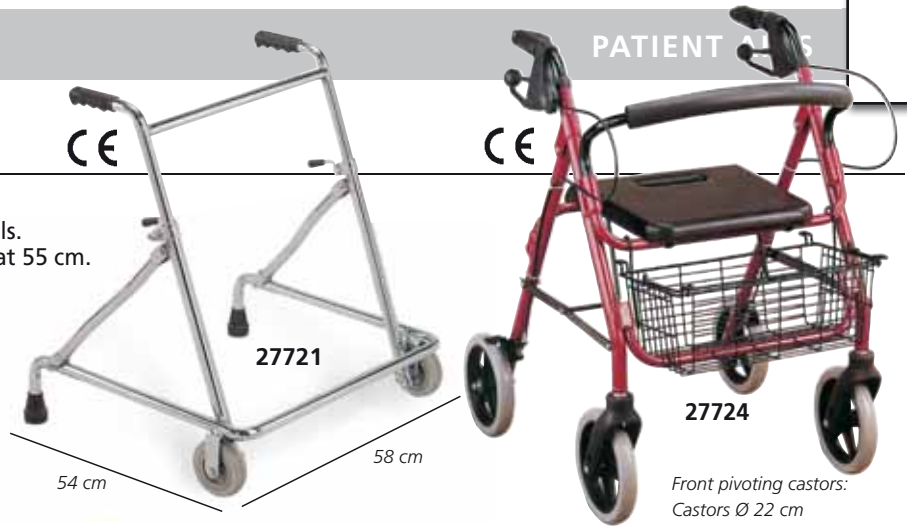

WALKER AID



• 27722 WALKER AID

Made of lightweight aluminium tubing, non-slip end and padded handgrips. Size 54x58 foldable with adjustable height 74.5-92 cm. Weight: 2.4 kg - Load: 100 kg.

• 27723 WALKER AID - WITH WHEEL

Adjustable height 85-100 cm. Size 49x44 cm.

• 27794 RUBBER TIP Ø 28 mm - grey for 27722-3 - bag of 5

PEDAL EXERCIZER

• 27738 PEDAL EXERCIZER

This pedal exerciser can be used on a table top as an arm exerciser or placed on the floor to exercise the legs. Hook-and-loop straps hold feet firmly to pedals. Convenient tension control knob allows a broad range of pedal resistance for a progressive exercise program.



WALKING STICKS-CRUTCHES



• 27780 ALUMINIUM WALKING STICK-silver

Adjustable height: 68-90.5 cm - weight 0.3 kg

• 27781 ALUMINIUM WALKING STICK - black

Adjustable height: 66-90.5 cm - weight 0.3 kg

• 27782 FOLDABLE ALUMINIUM WALKING STICK - coffee

Adjustable height: 82-92 cm - weight 0.3 kg

• 27790 QUADRIPODS - FOUR LEGS

Aluminium made, height adjustable 72.5-95 cm - 1.1 kg

• 27793 FOREARM CRUTCHES

Light weight aluminium tube. Anatomic handgrip.

Height adjustable 95.5-118.5 cm. Weight 0.95 kg. Max load 130 kg

• 27798 T-BAR CRUTCHES

Stainless steel type. Adjustable height: 118-130 cm. Max load 130 kg

• 27775 RACKER To reach and pick-up any object. Height: 82 cm

• 27795 RUBBER TIP Ø 19 mm - grey for 27780-2, 27793 - bag of 5

• 27796 RUBBER TIP Ø 16 mm - grey for 27790 - bag of 5

• 27797 RUBBER TIP Ø 22 mm - grey for 27798 - bag of 5

QUALITY STICKS MADE IN ITALY



WOOD STICKS:

• 27760 GIOTTO CURVED HANDLE - man - brown

• 27761 GIOTTO CURVED HANDLE - woman - brown

• 27762 RAFFAELLO "T" HANDLE - man - brown

• 27763 RAFFAELLO "T" HANDLE - woman - brown

• 27764 TIZIANO CURVED HANDLE - unisex - brown

• 27765 TIZIANO "T" HANDLE - unisex - brown

• 27768 LEONARDO METAL CURVED HANDLE - unisex - black

• 27769 LEONARDO METAL "T" HANDLE - unisex - black

• 27772 TIEPOLO CHROMED "DUCK" HANDLE - unisex - black

ALUMINIUM WHEELCHAIRS



CE

27719

• **27719 CLASSIC ALUMINIUM WHEELCHAIR**
Aluminium frame, flip-up desk armrest, detachable footrest, solid castor, quick release pneumatic rear wheel.



CE

27720

• **27720 KING ALUMINIUM WHEELCHAIR**
Adjustable aluminium chair frame, flip-up and height-adjustable armrest, detachable footrest. Quick release pneumatic rear wheels, solid castor, with carrying and anti wheel.



CE

LIGHT - ONLY 13 KG



27718

• **27718 LIGHT ALUMINIUM WHEELCHAIR**
Aluminium type with detachable armrest, footrest and brake folding, pneumatic rear wheel, solid castors.

POWERED WHEELCHAIR

• **27729 POWERED WHEELCHAIR**
Powered wheelchair, easy to use, for daily indoor use. The user can move the powered wheelchair forwards, backwards or turn round in the same place only by the controller. When the power is not enough, it can be used as a transit wheelchair. The seat cushion is detachable.
Main features:
Powder coated steel frame, adjustable armrest pad, detachable sideboards, detachable strong footrests, adjustable backrest inclination, rubber foam castors and drive wheels, with safety lamp.
Size: 1020x650x990 mm, weight: 50 kg (including batteries).
Load: 110 kg. Speed: 0-6 km/h.
Power: 2x12V=24V (20 km distance with fully charged battery)
Pivoting up to: 12°



COMMODOE CHAIRS

• **27701 COMMODOE CHAIR - stainless steel**
Stainless steel chair frame, with plastic detachable commode seat.
• **27702 COMMODOE WHEELCHAIR - painted**
Powder coated steel frame, detachable armrest and footrest, solid rear castors with lock, plastic commode seat.
• **27703 COMMODOE WHEELCHAIR - stainless steel**
Stainless steel chair frame, detachable armrest and footrest rear castors with lock, plastic commode seat.
• **27699 FOOTREST for 27702 - pair**
• **27700 FOOTREST for 27703 - pair**



27701



27702



27703

GIMA code	Chair	Chair	Wheelchair	Wheelchair	Wheelchair	Wheelchair	Wheelchair	Wheelchair	Wheelchair	Weight (Kg)	Dimensions (cm)	NW/GW
27702/3	55	44	12.5	12.5	51	92	90	44	42	100	54x53x53	17 19
27701	48	44			47	84		43	33	100	50x19x75	11

THERAPEUTIC ANTIDECUBITUS MATTRESSES UP TO LEVEL III

Code 28501

Antidecubitus therapy system for prevention and therapy of decubitus ulcers for low risk patients (Norton classification)



Made in Italy

- **28501 EURO BASIC ANTIDECUBITUS SYSTEM - 230V - 50 Hz**
Euro Basic is a 3 way alternating pressure therapy system suitable for low risk patients (Norton scale) with decubitus ulcers up to level II (NPUAP). Euro Basic is supplied with an electric compressor and an overlay mattress with 17 individually interchangeable inflatable elements made of PVC-PU; the head elements remain static to improve patients' comfort. The mattress can be easily connected to the compressor by the quick connectors and is quickly and safely fixed to the bed using the bed- straps with snap release buckles (supplied). The Euro Basic compressor incorporates a pressure regulator so that the correct inflation pressure can be set according to the patient's weight. It is designed to support a maximum of 150 Kg. The compressor has a visual alarm which indicates any system fault. The simple to operate CPR function ensures rapid deflation of the mattress for cardio pulmonary resuscitation. The mattress cover is breathable cover, waterproof and easily cleaned and can be machine washed up to 60°. System includes pump, overlay and cover. Available 110V - 60 Hz or 230V - 60 Hz on request

Code 28504

Prophylaxis and therapy of bedsores up to level II° (NPUAP classification). Medium risk patients (Norton classification)

Code 28506

Treatment and prevention of bedsores up to level III° (NPUAP classification). Medium-High risk patients (Norton classification)

Euro Soft Pro is an antidecubitus mattress replacement system with a digitally controlled compressor specifically designed to support patient weight up to 250 Kg.



Made in Italy

- **28504 EURO SOFT EVO ANTIDECUBITUS SYSTEM**
- **28506 EURO SOFT PRO ANTIDECUBITUS SYSTEM**

Mattress replacement system

- 20 individually interchangeable nylon-polyurethane inflatable elements, 18 cm deep
- Static head elements to improve patient's comfort
- Bi-elastic cover: breathable, bacteriostatic, antifungal, waterproof, fire proof, easy to clean, machine washable
- Polyurethane foam base, 5 cm deep, to ensure patient a sufficient support in case of emergencies
- CPR function for an easy and rapid deflation of the mattress (10 seconds)
- Transport facility: the mattress will remain inflated for 10-12 hours in case of power failure or during patient transport
- Inclination sensor to avoid sacral sinking
- Radiotranslucent surface enables examination without moving the patient
- Foam base fire proof (class 1IM) approved by the Ministry (approval number MI2435D20D1IM00002).

Cover

- Polyester fibre coated with PU. Elastic, water proof (ISO 811:1981), breathable, bacteriostatic, antifungal (ATCC 6538/ATCC 9642), easy to clean, machine washable up to 60°, fire proof (BS 7175:1989 - crib 5)

Compressor

- 3 way alternating pressure therapy, with static/dynamic function
- Digital control of mattress inflation pressure
- Max. inflation to assist patient ingress/egress and nursing
- Hour counter with alarm, indicating when maintenance is required.
- Audio/visual alarm with mute button, indicating any system fault
- Automatic key pad lock, to avoid any unauthorised changes of settings
- Acoustically insulated to ensure almost silent operation
- Size 27x23x13 cm, Weight 2 kg.

System includes pump, mattress, cover and inclination sensor.

- **28508 SPARE COMPRESSOR for 28504 and 28506**
Easy to install. Replacement is suggested approx every 24 months.

TECHNICAL SPECIFICATIONS

	28501	28504	28506
Prevention efficacy (Norton scale)	low risk	medium risk	medium-high risk
Prevention efficacy (NPUAP scale)	I level	II level	III level
Type of mattress	overlay mat	replacement mat	replacement mat
Mattress size - cm	200x90x13	200x90x18	200x90x23
Mattress weight - kg	7	7	9
Mattress material	nylon - PU	nylon - PU	nylon - PU
Max patient weight - kg	150	200	250
Deflation time (CPR)	10 sec.	10 sec.	10 sec.
Inflation time	25 min.	20 min.	30 min.
Cover material	nylon	polyester fibre with PU	
Power supply	220/230 V - 50 Hz		
Norms	CEI EN 60601-1:1998, CEI EN 60601-1-2:2001		



ANTI-DECUBITUS AIR MATTRESSES TO PREVENT SORES

- 28482 MASTER PUMP SYSTEM - with pressure adjustment
- 28483 KIT MASTER (AIR MATTRESS 28492 + MASTER PUMP 28482)

Made of stiff and strong ABS, it may be placed on the ground, as it has rubber feet, or hung at the edge of the bed. A control handle makes pressure regulation easier: air amount is adjustable to suit patient weight from 40 to 120 kg. The mattress should inflate totally in about 20 minutes. Air output 2.5 l/min. Master pump can be used with bubble or interchangeable section mattress. Flame retardant PVC.

Master pressure pump may be placed on the ground or hung at the edge of the bed.



Low noise, long lasting pump made with high quality materials

28482

TECHNICAL SPECIFICATIONS FOR MASTER PUMP - RINT PLUS PUMP

Size: 330x150x120 mm	Suitable for mattress: bubble or section
Weight: 1.95 Kg	Regulations: IEC-601-1 class II Type B
Operating voltage: 230 V - 50 Hz	Made in Italy
on request 110 V - 60 Hz	User manual in English Spanish, French, Italian, German
Power Consumption: 8 W - CLASS II	
Working pressure: 60-135 mm/Hg	

Patient Weight	Pressure
50kg - 110 lbs	40 mm/Hg
50-80 kg - 110-176 lbs	70 mm/Hg
80-150 kg - 176-330 lbs	100 mm/Hg



28492

MADE IN ITALY



28493

• 28492 AIR MATTRESS WITH BUBBLES

This air mattress is divided in bubble sections which inflate and deflate alternatively latitudinally through the alternative compressor. The continual movement given by the pressure pumped inside is useful to prevent minor bedsores problems.

• 28493 AIR MATTRESS - with interchangeable sections

It has 20 sections which can be inflated or deflated alternatively. Easy and practical to use, it is ideal for patients with minor bedsores problems.

Elements are jointed in-between and fixed to a PVC sheet.

• 28494 INTERCHANGEABLE CELL - spare for 28493

Press studs which attach each section to the supporting sheet allow an easy exchange in case of damage.

TECHNICAL SPECIFICATIONS MATTRESSES

	28487	28492	28493	28495
Total size:	190x85 cm	190x85 cm	190x85 cm	200x90x10.5 cm
Inflatable size:	160x85 cm	160x85 cm	160x85 cm	170x90 cm
Weight:	2.5 Kg	3 Kg	5.2 Kg	5.5 Kg
PVC thickness:	0.3 mm	0.4 mm	0.4 mm	0.4 mm
Number of elements:	-	-	20	17
Suitable for persons up to:	80/90 kg	80/90 kg	100 kg	120 kg

WELFARE KIT LAUNCH OFFER - MINIMUM 10 PCS

- 28489 - PUMP WITH PRESSURE ADJUSTMENT + BUBBLE MATTRESS (28487) Quality, reliable and economic kit with hook, made in Italy.
- 28487 WELFARE BUBBLE MATTRESS - beige

• 28495 RINT PLUS AIR MATTRESS - interchangeable cells

Rint Plus Mattress is made of a strong and soft PVC, realized with 17 big interchangeable bars. The integrated CPR function (Cardio Pulmonary Reanimation) allows the immediate emptying of the 17 bars. It is supplied with bi-elastic, liquid proof and permeable to water-steam cover. Made in Italy.

• 28498 RINT PLUS PUMP SYSTEM - with pressure adjustment Same as Master pump but with special valve for use with Rint Plus Mattress.

• 28496 INTERCHANGEABLE CELL - spare for 28495

• 28497 COVER - spare for 28495

28495



Thanks to height of cells, Rint Plus can be used also with medium risk patients.

ANTI-DECUBITUS PVC WATER MATTRESS

• 28499 WATER MATTRESS

It is made in PVC with a stabilising inflatable external structure to prevent excessive water "wave motion" and of 3 independent units. The water mattress distributes the total pressure on the whole supporting surface, making this less than the average arterial pressure, avoiding the exposure of the skin to pressure risk. Easy to use, it is placed over the usual mattress and it is washable. It is sold with a manual pump, 2 disinfectants and an inflation tube. Size: 190x85 cm. Weight 4 kg. PVC thickness 0,4 mm

Each water section measures 60x70 cm



28499

PVC rigid edge inflatable with air

RELAXING ARMCHAIR WITH 1 OR 2 ENGINES





45150
or 45155

10 COLOURS
AVAILABLE



45161
or 45166

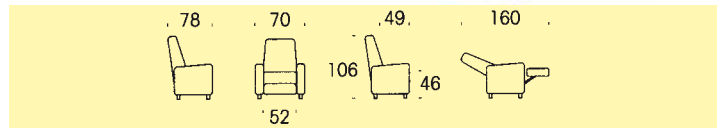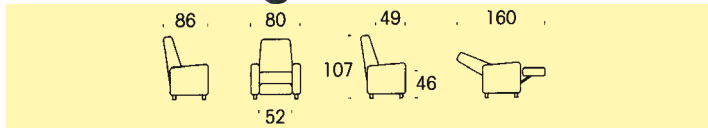

RELAXING ARMCHAIRS

Relaxing multifunctional armchair adaptable to any need. Available in 2 models (Beatrice and Betty), 1 or 2 engines and 10 colours. **One engine:** two button control for synchronized movement of leg section and back.

Two engines: four buttons control for independent movement of leg section and back.

Fabric 100% polyester. Load 120 Kg.

Lift motor 230V - 50Hz. Available on request 110V - 60 Hz or 230V - 60 Hz. Made in Italy.

GIMA code	BEATRICE ARMCHAIR
45150	Beatrice 1 engine - colour on request**
45151	Beatrice 1 engine - beige 2
45155	Beatrice 2 engines - colour on request**
45156	Beatrice 2 engines - beige 2

GIMA code	BETTY ARMCHAIR
45160	Betty 1 engine - colour on request**
45161	Betty 1 engine - blue 11
45165	Betty 2 engines - colour on request**
45166	Betty 2 engines - blue 11

** Product available in 15-20 days

HYDRAULIC AND ELECTRICAL LIFTERS

• **27749 HYDRAULIC LIFTER - max load 160 kg**
Lifting movement through a hydraulic pump operated by adjustable side control lever.
Descent movement by independent exhaust-valve.
Adjustable base 63 to 105 cm controlled by a manual device for an easy access to narrow places.
Painted steel frame.
Anti trauma push-handle.
2 rear castors Ø 100 mm, with independent brake and 2 front castors Ø 80 mm.
Equipped with universal sling.
Made in Italy.

• **27752 ELECTRICAL LIFTER 160 - max load 160 kg**
• **27753 ELECTRICAL LIFTER 200 - max load 200 kg**
Same as code 27749 but with electrical controls.
Up and down movements through an electric motor 24V D.C. and push button.
Battery fed with current transformer with an emergency push-button for the immediate stop of the movement and free up/down piston-stroke.

• **27757 UNIVERSAL SLING - spare for 27749, 27752**
Made of polypropylene, anti-allergic, non toxic, washable.

